


VascoPeel

VascoPeel is a radiopaque, pre-formable peel-away introducer set for use in the placement of: pacing leads, balloon catheters, temporary pacing catheters, drainage catheters, other catheters.

VascoPeel is designed to be removed without displacement of the introduced catheter. The risk of kinking or traumatic effects at the insertion point is minimized by its optimized design.

The material guarantees a good preformability and an excellent stability which is very important after removing the dilator out of the sheath.

A white screw fixation connects the dilator with the sheath and grants when advancing the introducer into the vessel that the dilator remains fixed with the sheath.

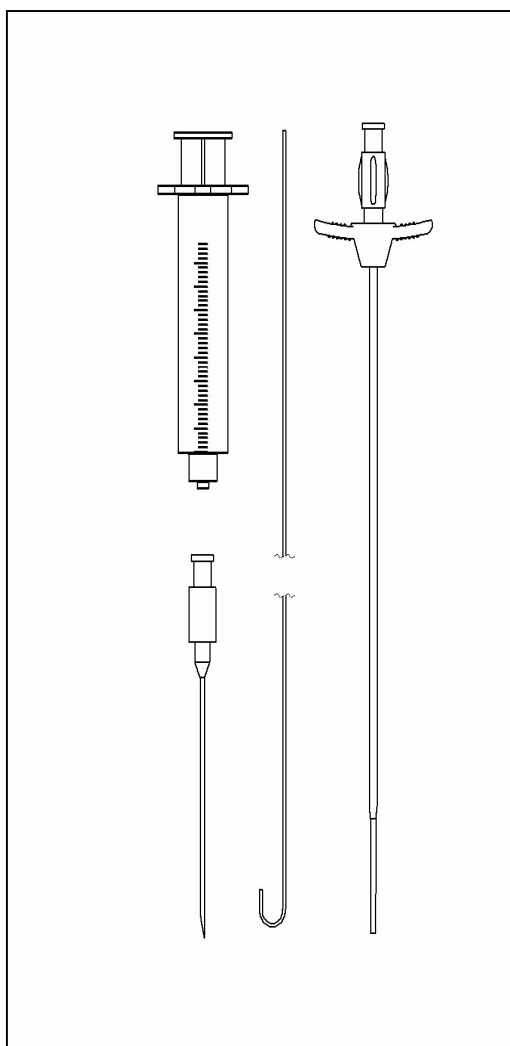

Technical Data:

Length (sheath):	14cm
Length (dilator):	17cm
French sizes:	5F to 12F

Package contents Set:

Set	Sheath/Dilator	Set mit 2 Sheath/Dilator
1 Peel-away Sheath introducer (14 cm)	1 Peel-away Sheath introducer (14 cm)	2 Peel-away Sheath introducer (14 cm)
1 Vessel dilator (17 cm)	1 Vessel dilator (17 cm)	2 Vessel dilator (17 cm)
1 0,035" x 45 cm safety J-tip guide wire		1 0,035" x 45 cm safety J-tip guide wire
1 thin-wall needle 18G		1 thin-wall needle 18G
1 LUER Slip syringe 10 ml		1 LUER Slip syringe 10 ml

Packaging Unit: 10 pcs. / Box

Artikel-Nummer:

Product	Set	Sheath/Dilator	Set with 2 Sheath/Dilator
VascoPeel 5F	211222	211431	-
VascoPeel 6F	211028	211105	-
VascoPeel 7F	211091	211106	211549
VascoPeel 8F	211029	211283	211094
VascoPeel 9F	211030	211198	211095
VascoPeel 10F	211036	211199	211096
VascoPeel 10,5F	211155	211075	211097
VascoPeel 11F	211035	211200	211098
VascoPeel 12F	211156	211201	211099

Accessories:	
J-Guidewire 0,035" x 45cm	208048

The products are double packaged and sterile (ETO-Sterilization).